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Multi-Quantum Beam High Voltage Electron Microscope Laboratory

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### The Laboratory of Multi-Quantum Beam High-Voltage Electron Microscope

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As a collaborative research and educational facility, we have developed a couple of high-voltage electron microscopes (hereinafter referred to as "HVEM") and various microscopes at Multi-Quantum Beam HVEM Laboratory of Hokkaido University. Since then, we have been supporting nano-analytical studies on surface and microstructure of materials.

In 1998 we introduced a multi-beam HVEM equipped with a double high-energy ion accelerator (ARM-1300) unlike any other that had ever existed in the world, then developed an advanced laser HVEM integrated with a pulsed-laser oscillator (Laser-HVEM) in 2007, and finally appended a laser irradiation device on a multi-beam HVEM to develop a multi-quantum beam HVEM in 2014. Therefore, we expect that our facilities will contribute to extremely unique and advanced research from the viewpoint of quantum beams of ion, electron and various optical quantum beams.

For the two research support projects, "The Advanced Research Infrastructure for Materials and Nanotechnology in Japan (ARIM Japan)" and "The Phase II of Inter-University Joint Research Centers of HVEM", we provide support services with the above mentioned equipment to hundreds of researchers in environment, energy and life-science relating to green-nano-life science every year.

In order to achieve our mission and further development, we still need broad public support. We hope that you take advantage of our world-leading facilities and research supporting systems, and greatly appreciate your further guidance and encouragement in the future.

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### A Common Facility for All Departments of the University

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The multi-beam HVEM is a system consisting of two ion accelerators (300kV and 400kV) and a 1300kV high resolution electron microscope with a point resolution of 0.117nm. Thus, we are now able to ionize and accelerate different ion species, to implant single or complex atoms simultaneously into materials and even to observe in situ their resolution at atomic level for the first time in the world.

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
Multi-Quantum Beam HVEM Laboratory of Hokkaido University participates in “The Advanced Research Infrastructure for Materials and Nanotechnology in Japan of the Ministry of Education, Culture, Sports, Science and Technology (MEXT)”, and we will improve the joint-use system of the facilities, promote mutual network systems and support nanotechnology research involving high-voltage electron microscopy.

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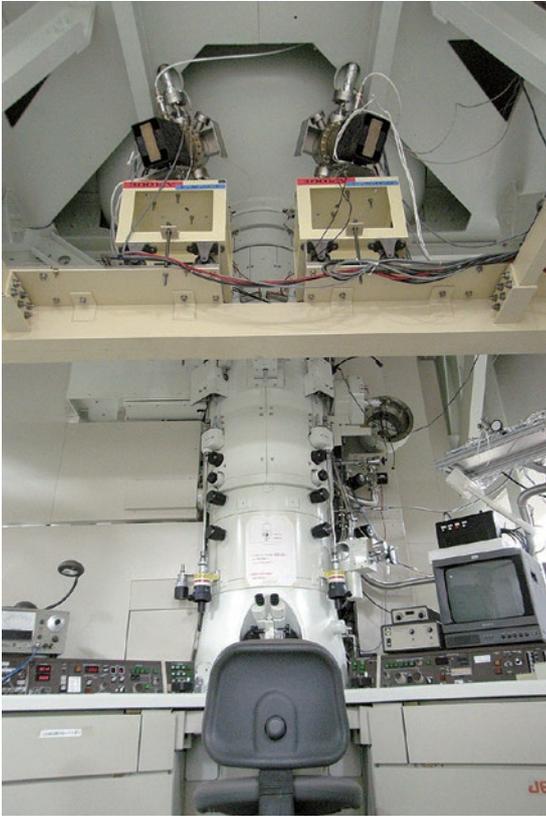### ▶▶主な用途 Primary use

超高分解能構造像・格子像の観察  
Observation of ultra-high resolution structure image and lattice image  
厚膜試料の観察  
Observation of thick film sample  
試料加熱・冷却過程その場観察  
In-situ observation of microstructural change during heating/  
cooling processes  
イオン・レーザー・電子線照射その場観察  
In-situ observation of ion implantation/radiation damage behaviors

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### ▶主な用途 Primary use

透過型電子顕微鏡観察用薄膜試料制作  
Preparation of thin film sample for transmission electron microscopy  
SIM像観察  
Observation of SIM image

## 電界放射型 200kV 透過型電子顕微鏡 JEM-2010F 200kV Field-Emission Transmission Electron Microscope

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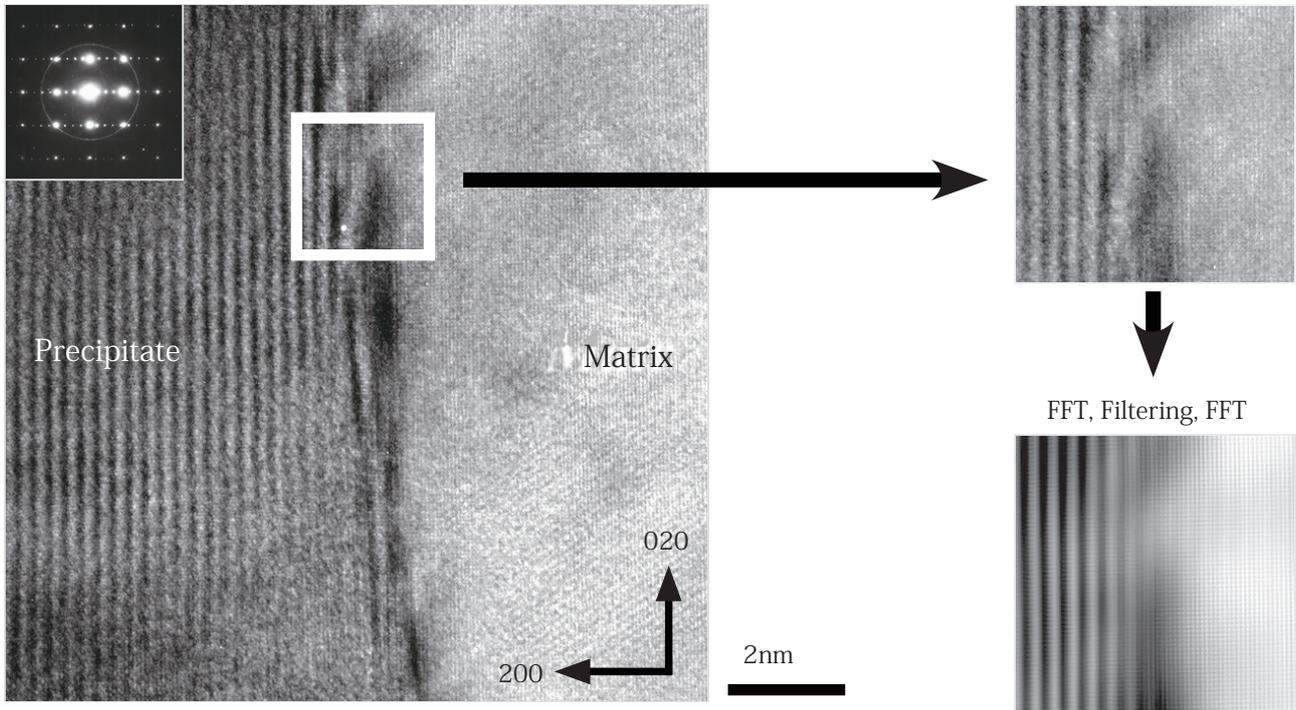
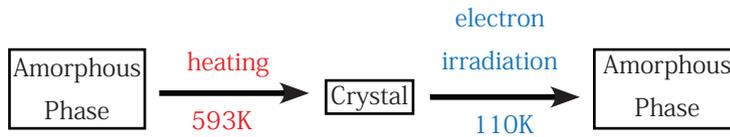
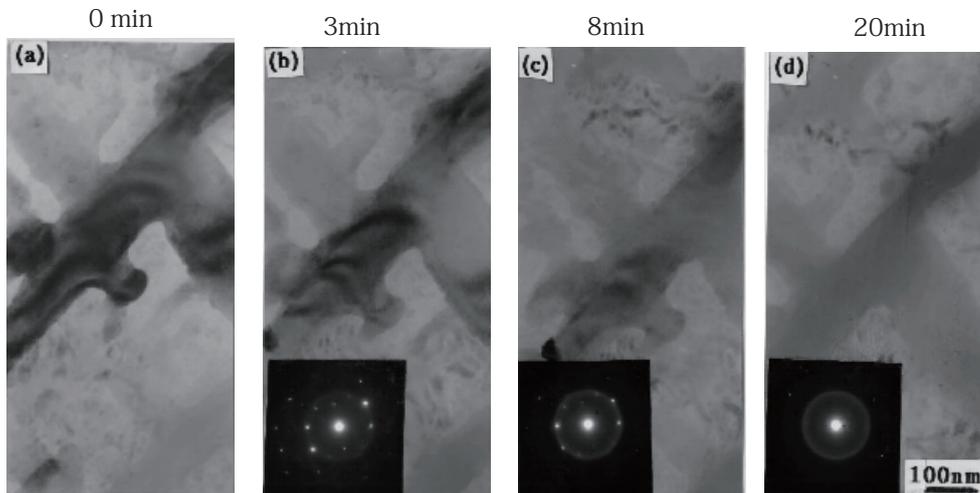
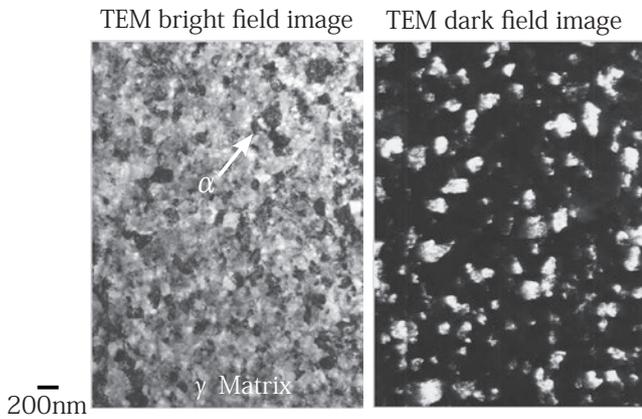
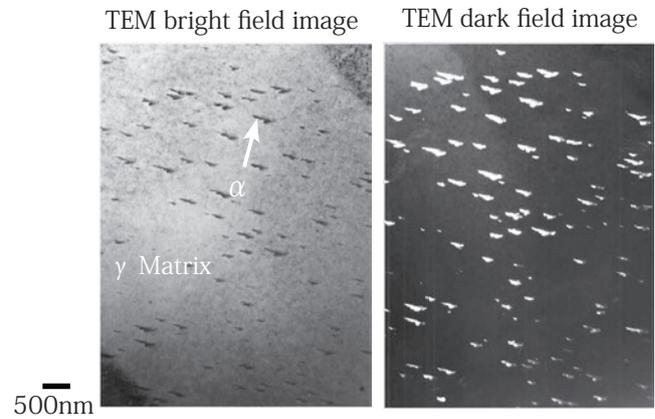
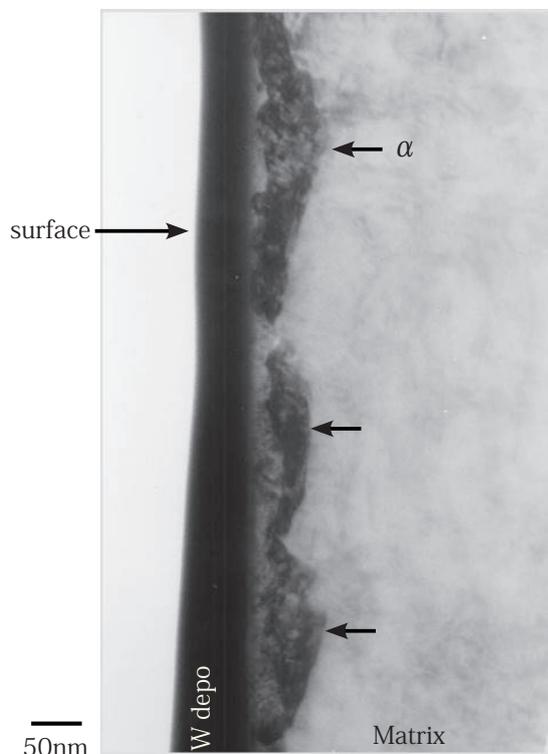
Cross-sectional TEM image of phase transformation

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
Now it is possible to manufacture materials with different structures of phases and atoms by using the ion accelerator equipped on the multi-beam HVEM, while it is not obtainable by existing techniques. For example, implanting various metallic elements to aluminum alloys by ion implantation disperses nonequilibrium phase deposition finely, and thus it can fortify only the ion implanted layer.


Also, it enables to form a hard martensitic phase on the stainless steel for the construction industry. Therefore, we expect that this new surface modification technology will create further applications in the future.

H. Kinoshita et al., Mater. Chem. Phys., 81 (2003)  
Dwi Gushiono et al., Mater. Trans., 45 (2004)

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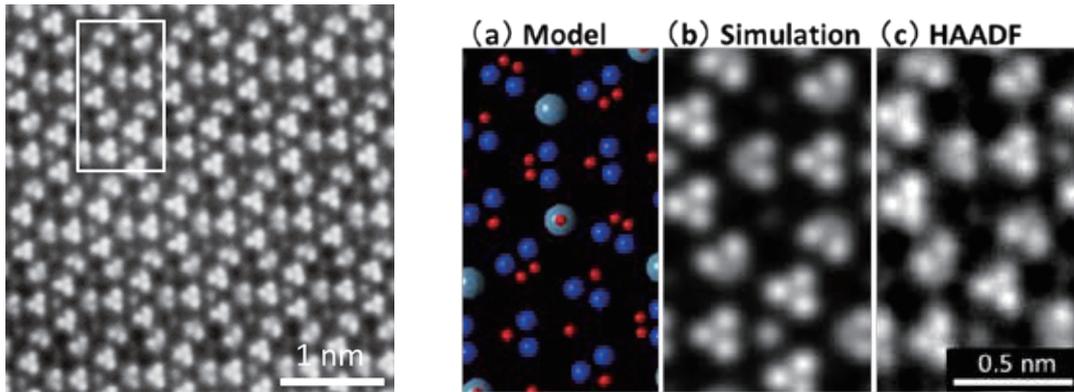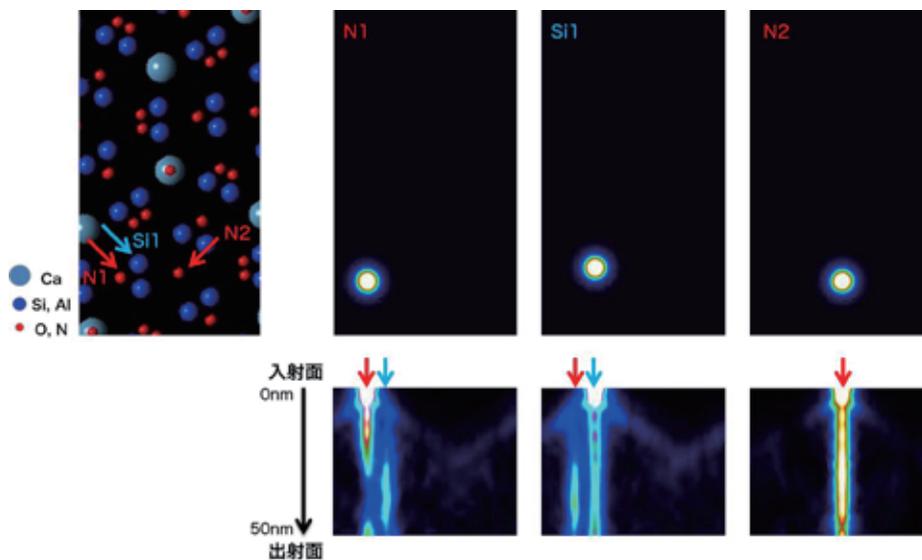
Due to the miniaturization of devices in recent years, it is essential to carry out research and development of advanced materials for their atomic structural evaluation and chemical bonding state analysis at sub-nanometer scale. For this purpose, we carry out analysis of advanced device materials by combining a multi-beam HVEM and various analytical electron microscopes.

For example we have observed directly atomic structures of metal/ceramic interface under multi-beam HVEM and also found out specific electronic states at the interface by a combination of electron energy-loss spectroscopy (EELS) and first-principle electronic structure calculations. In addition, we have clarified the dopant distribution of advanced ceramic materials by using the spherical aberration corrected STEM and the multi-slice simulation technique.

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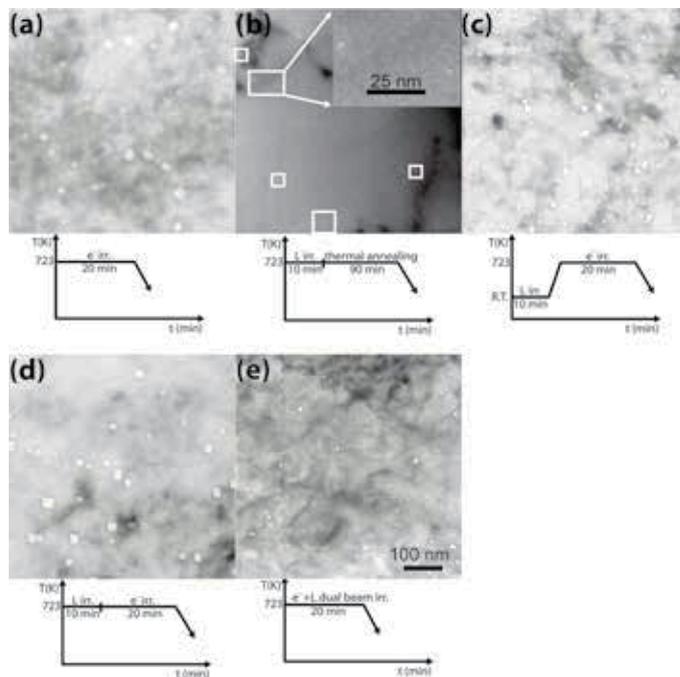
Void formation under laser-electron irradiation (a) electron irradiation (b) laser irradiation (c,d) laser irradiation followed by electron irradiation (e) simultaneous laser and electron dual-beam irradiation

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
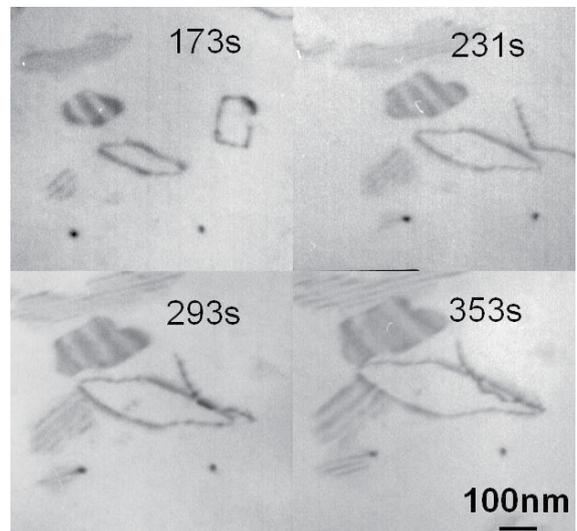
We developed a multi-beam HVEM connected with two ion accelerators in 1998. Subsequently, we developed another HVEM with an integrated short-pulsed laser in 2007, then developed a single ion laser radiation device equipped with an ion accelerator whose vacuum chamber was provided with a short-pulse laser unit in 2009. Currently we have accomplished to install a multi-beam HVEM applied to the three quantum beams (laser, electron and ion). When we operate it, all the representative elementary particles, namely photon (laser beam), lepton (electron) and hadron (ion) will become a complete set. Thus, we will be able to implement simulation experiments that has been unfeasible, such as dealing with damage materials exposed to high temperature plasma in a nuclear fusion reactor where gamma rays (electromagnetic radiation), beta rays and neutron/ion beams are coexisting and competing. Since it makes us possible to conduct a nano-region observation and an in situ observation simultaneously under electron microscope, we can expect further research promotion on high time and high spatial resolution analysis for green-nanotech-energies.

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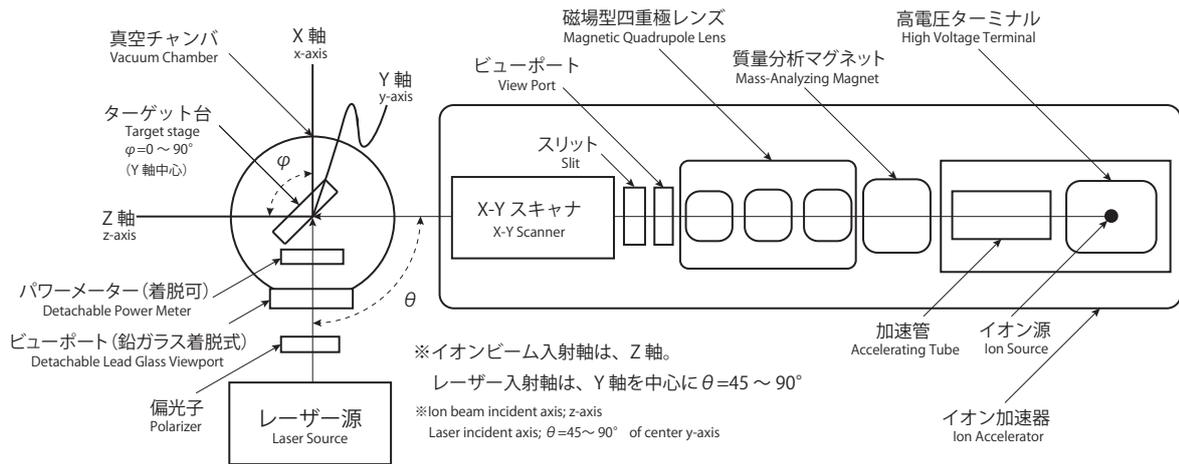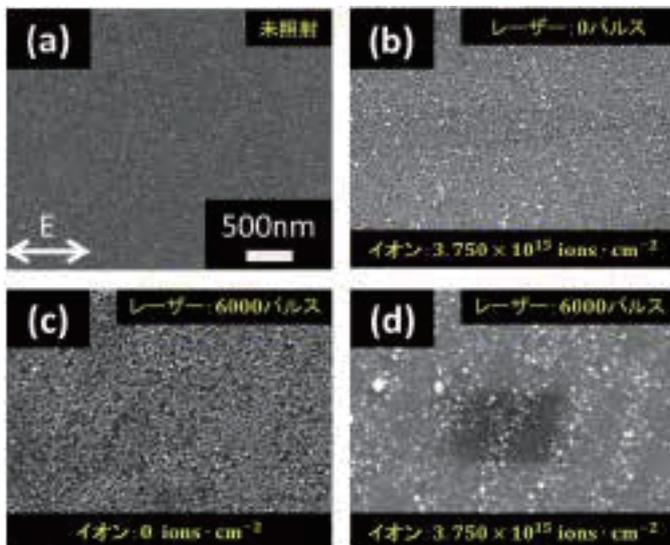
SEM images of before and after irradiation. (a) non-irradiated.  
(b) After ion irradiation. (c) After laser irradiation.  
(d) After ion and laser co-irradiation.

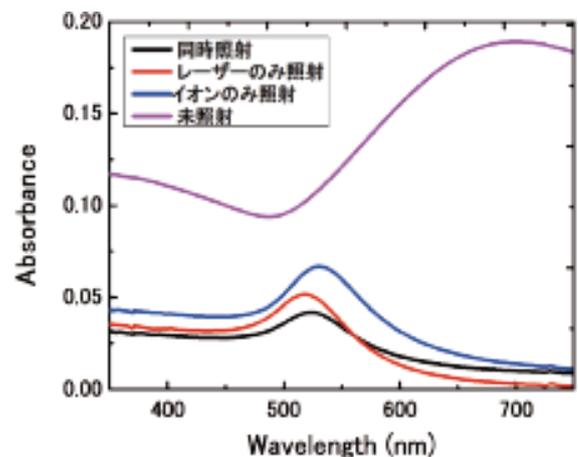照射領域の吸収スペクトル  
Absorption Spectrum of Irradiated Specimen

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It has become clear that simultaneous irradiation with ion and short-pulse laser can be applied to create nanoparticles of noble metals. Therefore, our HVEM laboratory proceeds to investigate the mechanism of new nanostructure formations by a combination of the ion accelerator and the short-pulse laser.

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

The wavelength of electrons is very short due to the acceleration of the multi-beam HVEM. Consequently it gives high resolution images and makes electron beams transmit easily through the sample. In particular, when we observe materials that are weak against irradiation of electron beams, such as biological materials at high resolution under a conventional electron microscope (accelerating voltage 200kV), the sample is damaged by electron beam heatings. On the contrary, we can observe such materials at atomic level using the multi-beam HVEM without heating.

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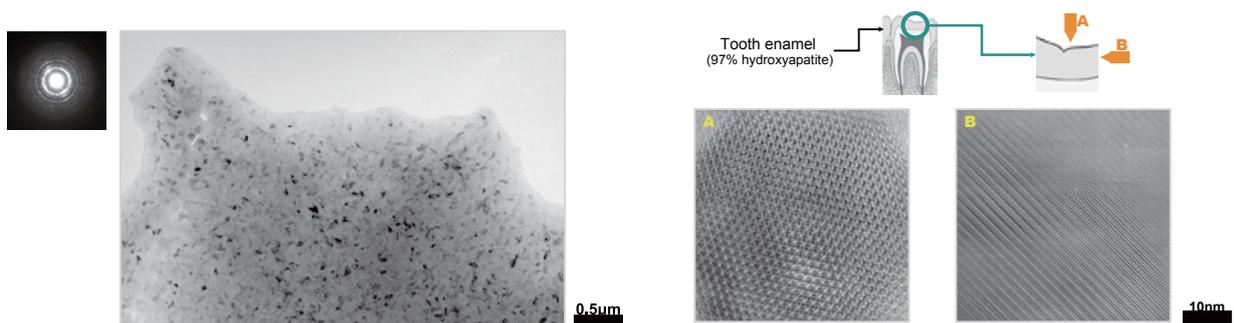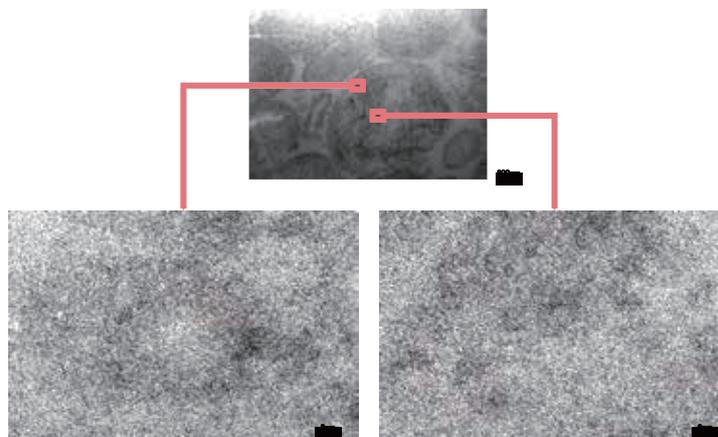Any personnel intending to use the equipment in our facilities for ARIM Japan Project must follow the following procedures:

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Anyone intending to use the equipment in our facilities must follow the following procedures:

1. 利用申請書の提出 (はじめて当施設をご利用の方) / Submission of the form for use of the facilities (for first-time users)  
利用希望者は利用申請書に必要事項を記入の上、複合量子ビーム超高压顕微解析研究室へご提出ください。  
Please submit the application form by e-mail to Multi-Quantum Beam HVEM Laboratory.
2. 設備利用予約 / Booking for the use of the facilities  
利用申請後、電話またはメールにて、希望利用設備と日時をお知らせください。  
Please inform the person-in charge via phone or e-mail about your booking ( which facility to use, date and time) after the submission of the form.

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Usage fees will be billed through the Open Facility service. Please apply to use the Open Facility service from the website below.  
<https://www.gfc.hokudai.ac.jp/>

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# User's Guide and Access to the Facilities

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### ▶ 空港～JR 利用 / From New Chitose Airport by JR trains

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Travel to Sapporo Station by Rapid Airport train. Then, travel to Kita 12 Jo or Kita 18 Jo Subway Station of the Namboku line.

It is a tenminute walk to the Laboratory of HVEM.  
(Please see the map.)

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Buses from New Chitose Airport/


By Hokuto Kotsu Bus: Get off at the South Gate of Sapporo Station.

By Chuo Bus : Get off at Keio Plaza Hotel.

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After getting off the bus, take a subway of the Namboku line from Sapporo Subway Station to Kita 12 Jo or Kita 18 Jo. It is a ten-minute walk to the Laboratory of HVEM.

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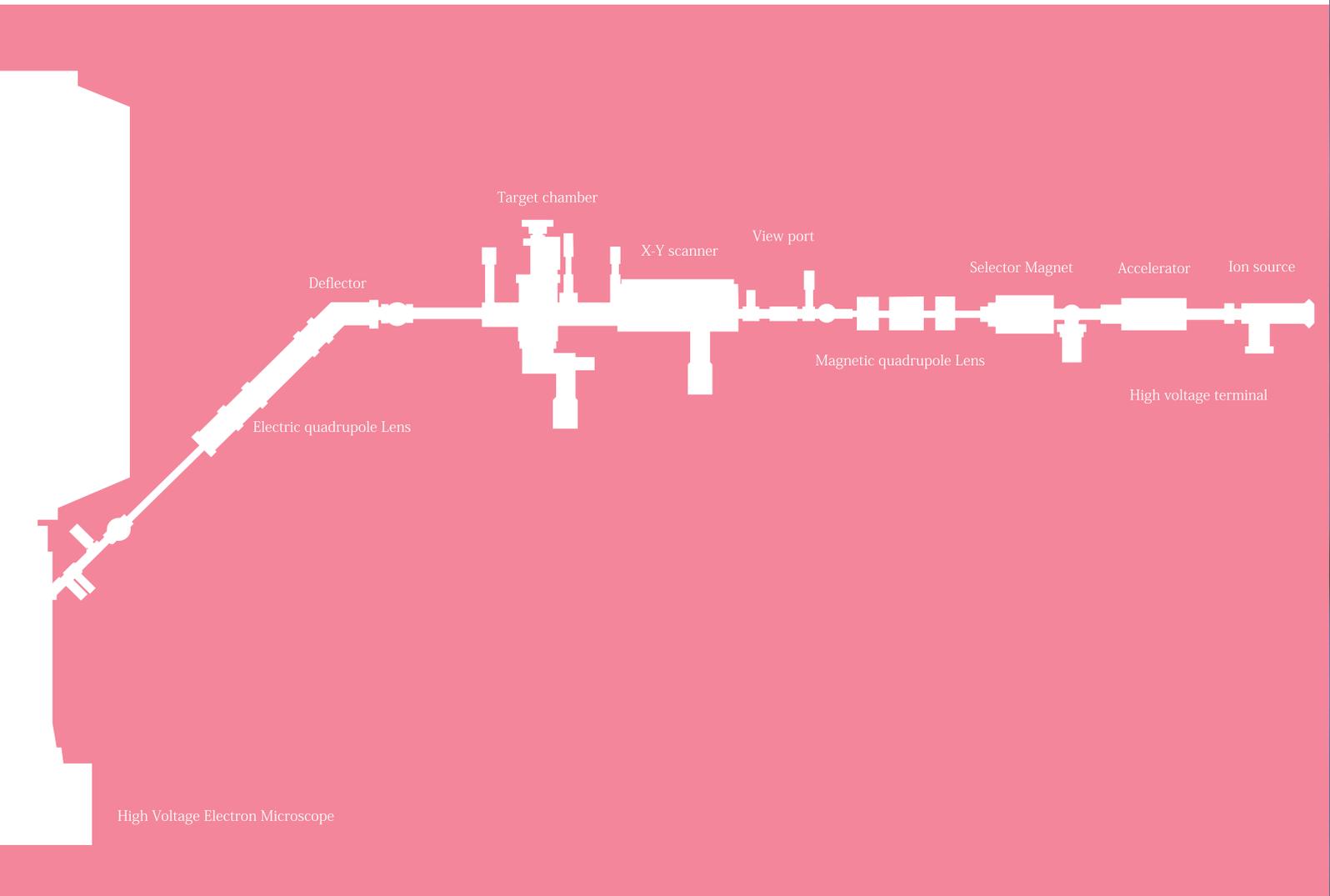Multi-Quantum Beam HVEM Laboratory

Center for Advanced Research of Energy and Materials, Faculty of Engineering of Hokkaido University

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